


ANALYTICAL REPORT

4F-CUMYL-5F-PINACA (C₂₂H₂₅F₂N₃O)

1-(5-fluoropentyl)-N-[2-(4-fluorophenyl)propan-2-yl]-1H-indazole-3-carboxamide

Remark – other active cpd. detected **none**

| | |
|-----------------------------|-----------------------------|
| Sample ID: | 1946-18 |
| Sample description: | powder - granulated - white |
| Sample type: | RM-reference material |
| Comments: | Chiron AS Lot20662, |
| Date of entry (DD/MM/YYYY): | 01/06/2018 |
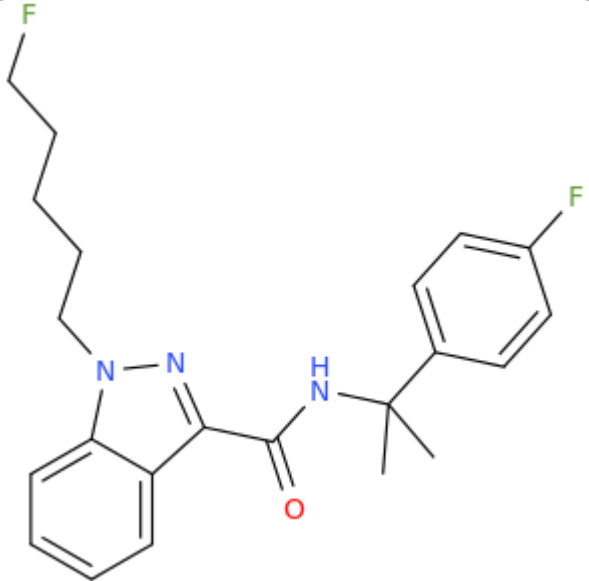
| | |
|---|---|
| Substance identified-structure ¹ (base form) |  |
| Systematic name: | 1-(5-fluoropentyl)-N-[2-(4-fluorophenyl)propan-2-yl]-1H-indazole-3-carboxamide |
| Other names: | SGT-65 |
| Formula (per base form) | C ₂₂ H ₂₅ F ₂ N ₃ O |
| M _w (g/mol) | 385,46 |
| Salt form: | base |
| StdInChIKey (per base form) | WAXVVLQFKMMDHR-UHFFFAOYSA-N |
| Other active cpd. detected | none |
| Add.info (purity..) | pure by GC-MS, HPLC-TOF |

¹ Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

| date | comments (explanation) |
|------|------------------------|
| | |
| | |
| | |
| | |

Supporting information

| Analytical technique: | applied | remarks |
|-------------------------|---------|---|
| GC-MS (EI ionization) | + | NFL GC-RT (min): 11,76 BP(1): 233; BP(2): 145,BP(3) :370, |
| FTIR-ATR | + | direct measurement |
| GC-IR (condensed phase) | + | always as base form |
| HPLC-TOF | + | exact mass theoretical: 385,1966 / measured Δ ppm: -1,61 |

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 μ m. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 7.1 min, then heating at 50 °C/min up to 325 °C and finally 6.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm^{-1} ; resolution 4 cm^{-1}

3. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above (1). Split MS: IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm^{-1} .

4. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; injection volume 1 μ l. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N₂) and sheath temperature 325 °C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.


ANALYTICAL RESULTS

MS (EI)

Abundance

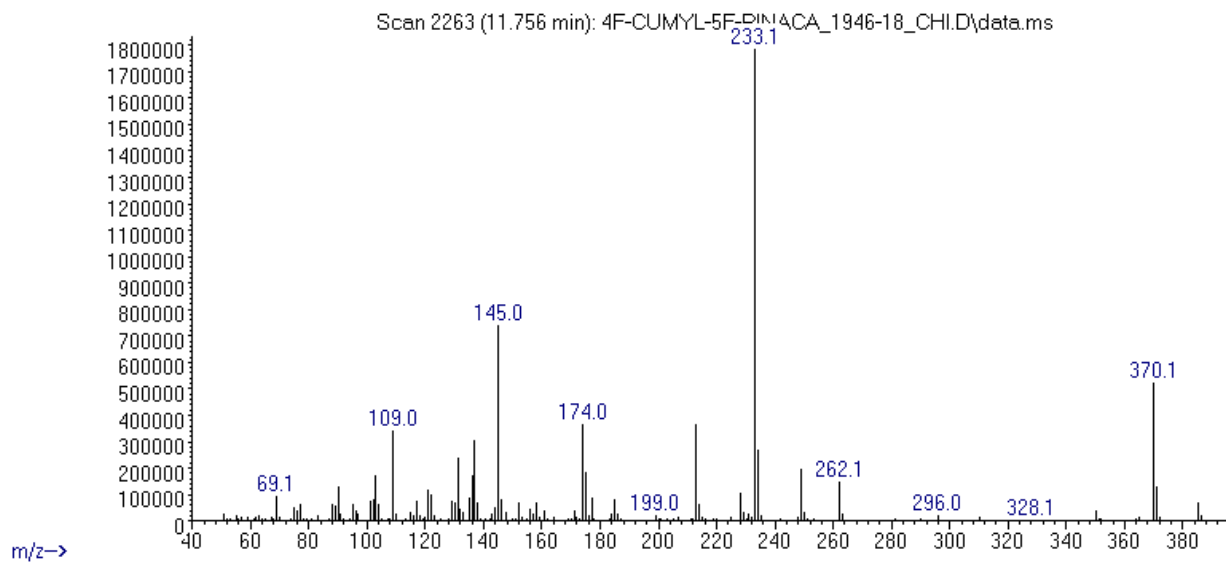

FTIR-ATR

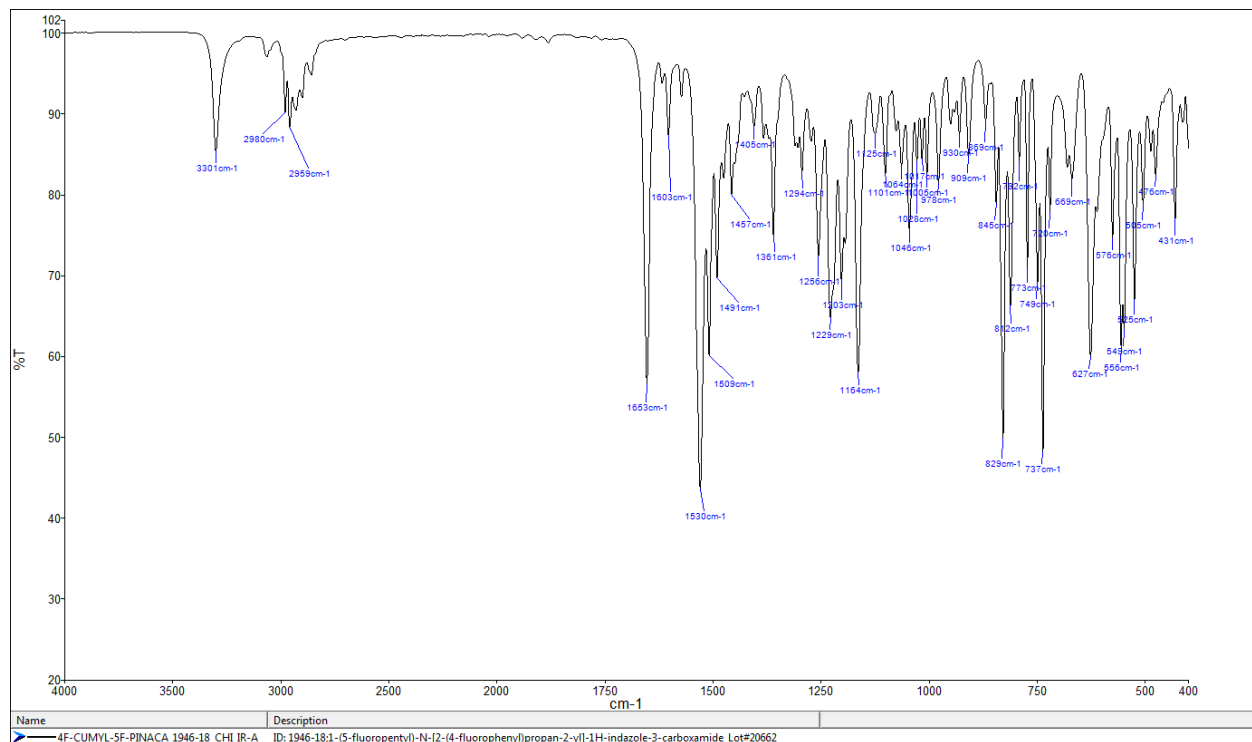

IR (condensed phase – after chromatographic separation)

