


ANALYTICAL REPORT

9(R)-HHC-P acetate (C₂₅H₃₈O₃)

(6aR,9R,10aR)-3-heptyl-6,6,9-trimethyl-6H,6aH,7H,8H,9H,10H,10aH-benzo[c]isochromen-1-yl acetate

Remark – other active cpd. detected: none

Sample ID:	3242-23
Sample description:	liquid - transparent
Sample type:	RM-reference material
Comments:	CAY Lot#0663993-4,
Date of entry (DD/MM/YYYY):	26/07/2023
web link	http://www.policija.si/apps/nfl_response_web/seznam.php

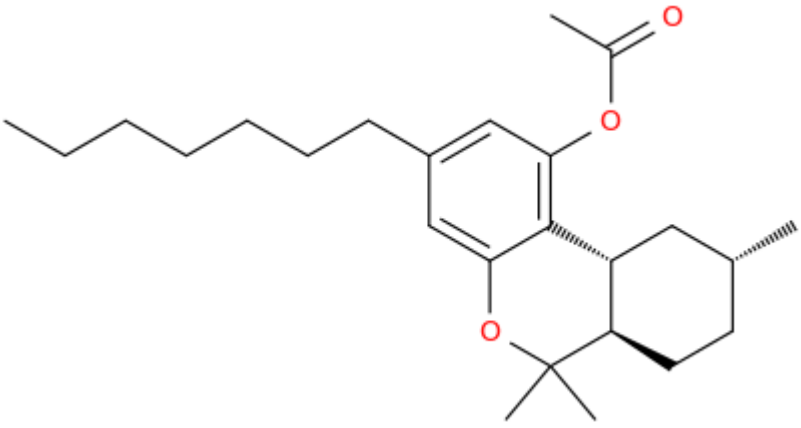Substance identified-structure ¹ (base form)	
Systematic name:	(6aR,9R,10aR)-3-heptyl-6,6,9-trimethyl-6H,6aH,7H,8H,9H,10H,10aH-benzo[c]isochromen-1-yl acetate
Other names:	9(R)-HHCP Acetate; 9(R)-Hexahydrocannabiphorol acetate; (6aR,9R,10aR)-3-heptyl-6,6,9-trimethyl-6a,7,8,9,10,10a-hexahydro-6H-benzo[c]chromen-1-yl, acetate
Formula (per base form)	C ₂₅ H ₃₈ O ₃
M _w (g/mol)	386,58
Salt form:	base
StdInChIKey (per base form)	LRVWJGANNJLTC-DUXKKGJEZSA-N
Other active cpd. detected	none
Add.info (purity..)	≥98%; 5mg/0,5 ml acetonitrile solution

¹ Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 10,19 BP(1): 260; BP(2): 301,BP(3) :344,
FTIR-ATR	-	direct measurement
GC-IR (solid phase)	+	always as base form
HPLC-TOF	+	exact mass theoretical: 386,2821 / measured Δ ppm: -2,53

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 μ l and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 μ m. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 9.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm^{-1} ; resolution 4 cm^{-1}

3. GC- (MS)-IR solid phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 μ l and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above (1). Split MS: IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm^{-1} .

4. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; injection volume 1 μ l. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N₂) and sheath temperature 325 °C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.

ANALYTICAL RESULTS

MS (EI)

Abundance

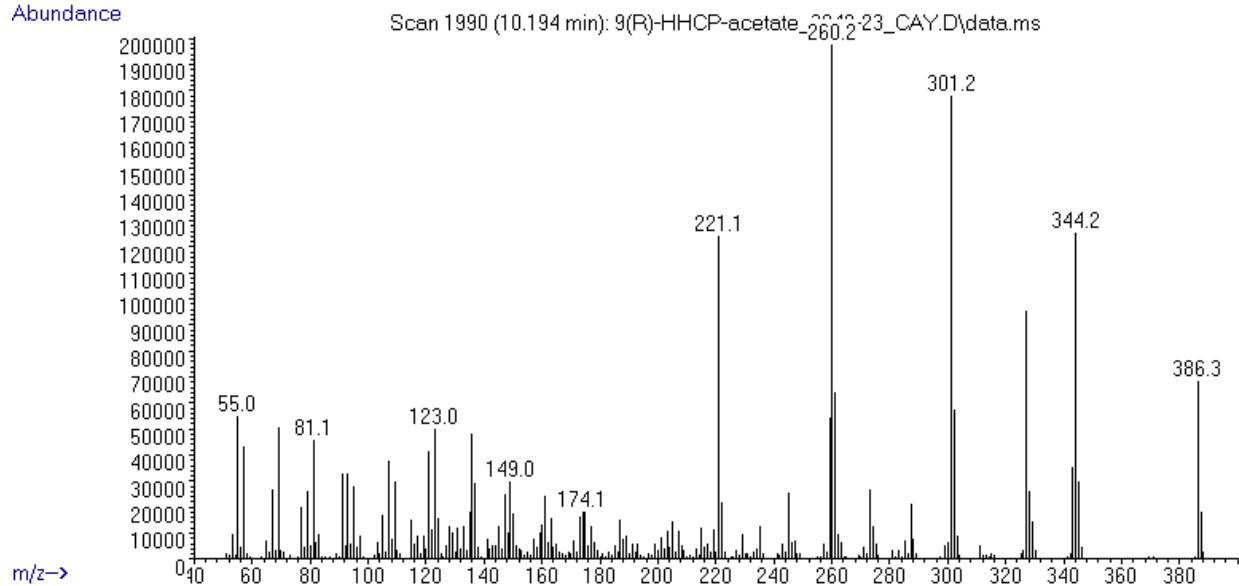

IR- (solid phase – after chromatographic separation) - spectrum per base form

