


ANALYTICAL REPORT

Clonitazene (C₂₀H₂₃ClN₄O₂)

(2-[2-[(4-chlorophenyl)methyl]-5-nitro-1H-1,3-benzodiazol-1-yl]ethyl)diethylamine

Remark – other active cpd. detected: none

Sample ID:	3044-21
Sample description:	powder - white
Sample type:	RM-reference material
Comments:	Chiron Lot#26173,
Date of entry (DD/MM/YYYY):	10/01/2024
web link	http://www.policija.si/apps/nfl_response_web/seznam.php

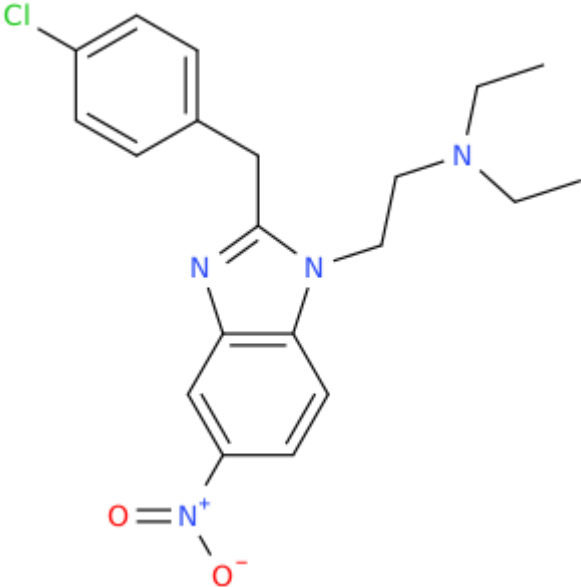Substance identified-structure ¹ (base form)	
Systematic name:	(2-[2-[(4-chlorophenyl)methyl]-5-nitro-1H-1,3-benzodiazol-1-yl]ethyl)diethylamine
Other names:	2-[2-[(4-chlorophenyl)methyl]-5-nitrobenzimidazol-1-yl]-N,N-diethylethanamine; 2-(p-Chlorobenzyl)-1-(2-diethylaminoethyl)-5-nitrobenzimidazole; Clobedolum
Formula (per base form)	C ₂₀ H ₂₃ ClN ₄ O ₂
M _w (g/mol)	386,88
Salt form:	HCl
StdInChIKey (per base form)	GPZLDQAEBHMPG-UHFFFAOYSA-N
Other active cpd. detected	none

¹ Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Add.info (purity..)	99,9%
---------------------	-------

Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 16,1 BP(1): 86; BP(2): 125,BP(3) :87,
FTIR-ATR	+	direct measurement
GC-IR (solid phase)	+	always as base form
HPLC-TOF	+	exact mass theoretical: 386,151 / measured Δ ppm: -0,27

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 μ l and split mode (1:50). Injector temperature: 280 $^{\circ}$ C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 μ m. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 $^{\circ}$ C for 1 min, followed by heating up to 190 $^{\circ}$ C at rate 8 $^{\circ}$ C/min, then heating up to 293 $^{\circ}$ C at a rate of 18 $^{\circ}$ C/min, hold for 6.1 min, then heating at 50 $^{\circ}$ C/min up to 325 $^{\circ}$ C and finally 9.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235 $^{\circ}$ C, source and quadropole temperatures 280 $^{\circ}$ C and 180 $^{\circ}$ C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm^{-1} ; resolution 4 cm^{-1}

3. GC- (MS)-IR solid phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 μ l and split mode (1:5). Injector temperature 280 $^{\circ}$ C. Chromatographic separation as above **(1)**. Split MS: IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235 $^{\circ}$ C, source and quadropole temperatures 280 $^{\circ}$ C and 180 $^{\circ}$ C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm^{-1} .

4. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; Injection volume 1 μ l. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N₂) and sheath temperature 325 $^{\circ}$ C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.


ANALYTICAL RESULTS

MS (EI)

Abundance

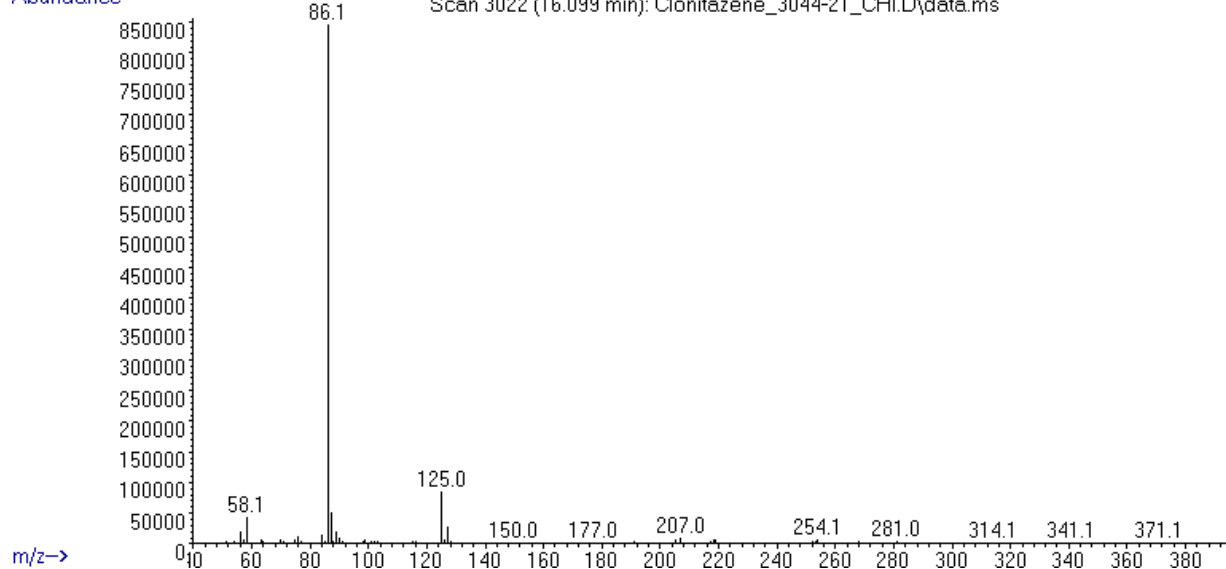

FTIR-ATR (direct measurement – sample as received)

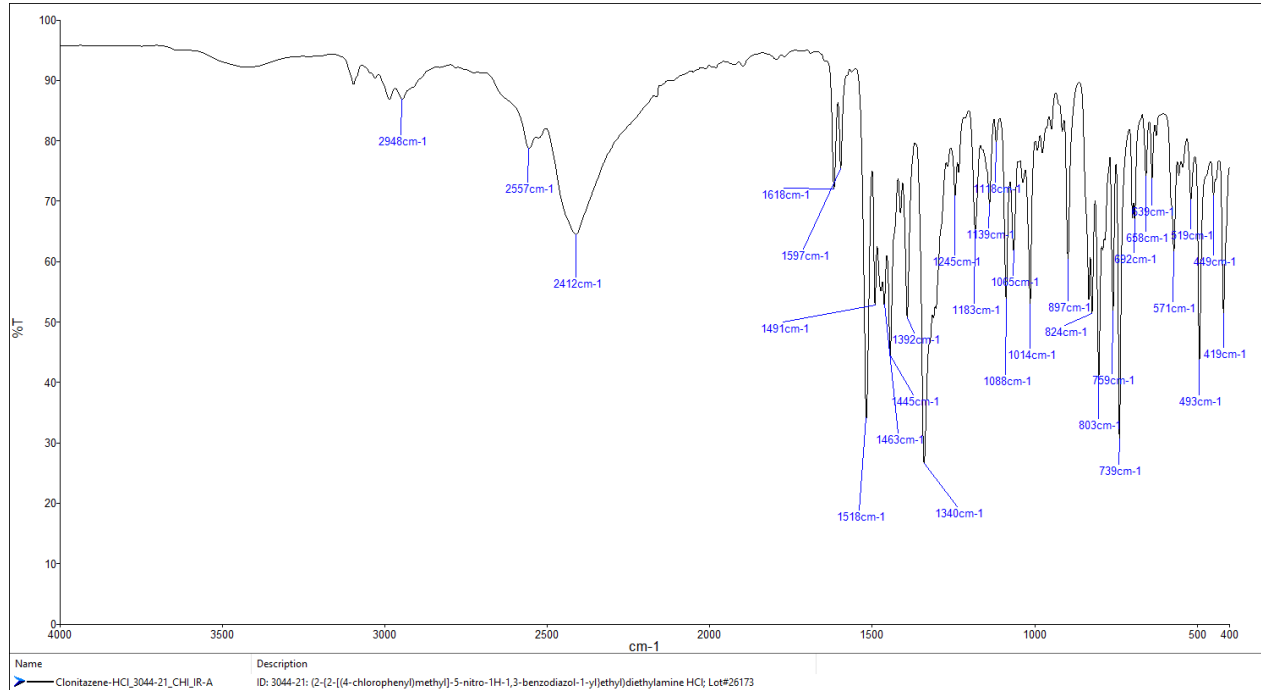

IR- (solid phase – after chromatographic separation) - spectrum per base form

