


ANALYTICAL REPORT

N-Methylethylone (C₁₃H₁₇NO₃)

1-(2H-1,3-benzodioxol-5-yl)-2-[ethyl(methyl)amino]propan-1-one

Remark – other active cpd. detected: **none**

Sample ID:	1919-18
Sample description:	powder - white off
Sample type:	RM-reference material
Comments:	CAY Lot#0519716-4,
Date of entry (DD/MM/YYYY):	09/04/2018

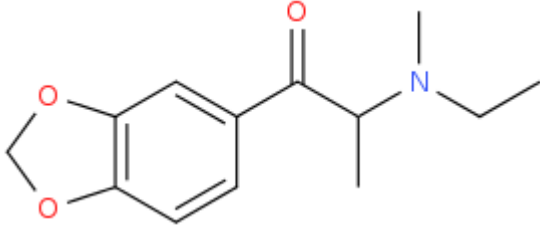Substance identified- structure ¹ (base form)			
Systematic name:	1-(2H-1,3-benzodioxol-5-yl)-2-[ethyl(methyl)amino]propan-1-one		
Other names:	N-Ethylmethylone;	3,4-Methylenedioxy-N-ethyl-N-methylcathinone;	1-
	(benzo[d][1,3]dioxol-5-yl)-2-(ethyl(methyl)amino)propan-1-one; MDMEC		
Formula (per base form)	C ₁₃ H ₁₇ NO ₃		
M _w (g/mol)	235,28		
Salt form:	HCl		
StdInChIKey (per base form)	AKDKKEIHDXROHU-UHFFFAOYSA-N		
Other active cpd. detected	none		
Add.info (purity..)	100 %		

¹ Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 5,89 BP(1): 86; BP(2): 58,BP(3) :87,
FTIR-ATR	+	direct measurement
GC-IR (condensed phase)	+	always as base form

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 µm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 7.1 min, then heating at 50 °C/min up to 325 °C and finally 6.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

3. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above **(1)**. Split MS : IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm⁻¹.

4. HPLC-TOF for exact monoisotopic mass and empirical formula control - results are not shown in the report.


ANALYTICAL RESULTS

MS (EI)

Abundance

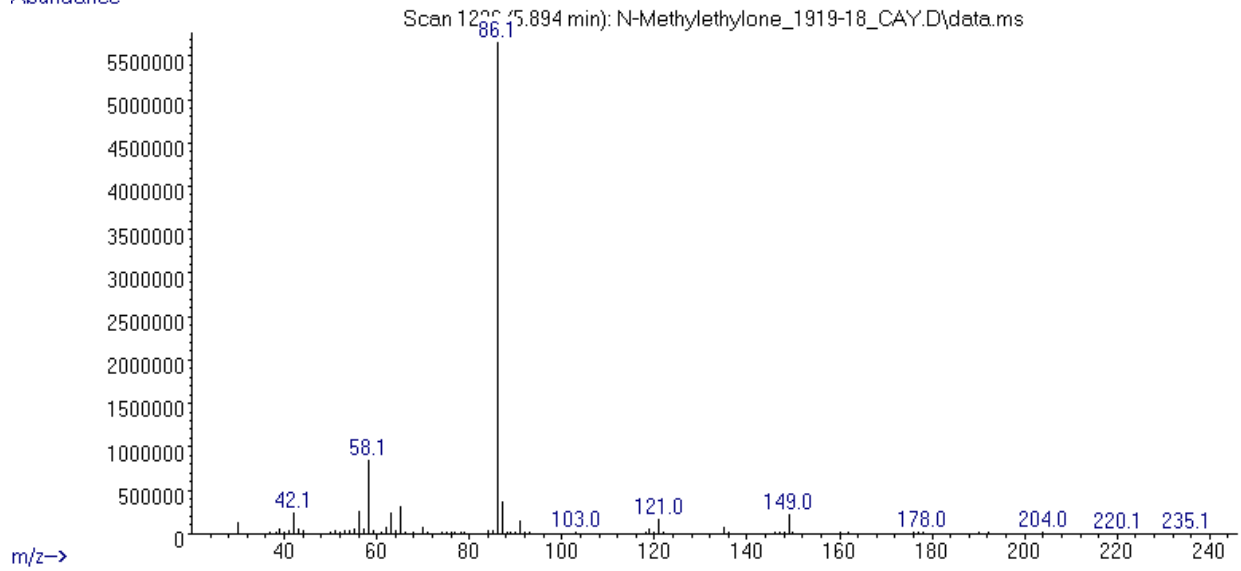

FTIR-ATR - sample as received

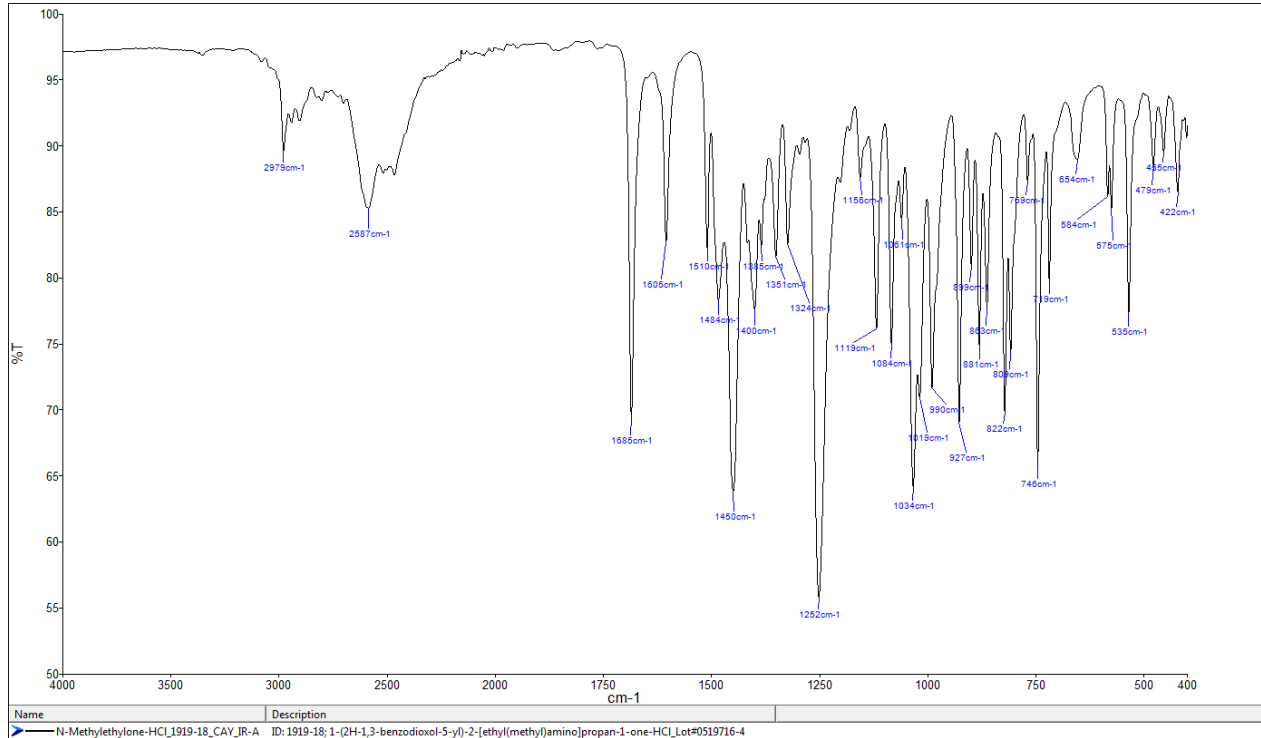

IR (condensed phase – after chromatographic separation)

